

2nd Generation



NObreath[®] V2

- ✓ **FDA Cleared**
- ✓ **Conforms to ERS/ATS Guidelines**
- ✓ **Recommended by NICE**



호산구성 기도염증이 진행 중인 천식 및 기도 감염에서
호기산화질소의 수치가 증가합니다.

NObreath[®]V2는 아토피성 천식 및 호산구성 염증의
증가로 인한 천식의 진단에 유용합니다.

DNAMED 동남메디칼

경기도 고양시 일산동구 대산로 18 서교빌딩 2층
Tel: 031) 908-7588 <http://www.dnamed.co.kr>



Measure breath nitric oxide for airway inflammation with the NObreath® FeNO monitor

NObreath® V2



Clinically proven

FDA 승인
ERS/ATS 지침 준수
NICE 추천




FeNO testing without limits

적절한 사용법과
유지보수로 무제한 검사



Easy to use

화면을 보고 따라하면 검사완료

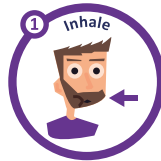

Adult & child testing



Technical specifications

Concentration range		5-500ppb
Display		Full colour touchscreen
Detection principle		Electrochemical sensor
Accuracy		±5ppb of measured value ≤ 50ppb ±10% of measured value > 50ppb
Breath test time		Adult: 12 seconds Child: 10 seconds Ambient: 30 seconds
Temperature	Operating	15-30°C
	Storage/transport	0-50°C
Dimensions		Approx. 90 x 159 x 59 mm
Weight		Approx. 400g

수입원



경기도 고양시 일산동구 대산로 18 서교빌딩 2층
Tel: 031) 908-7588 / Fax: 031) 908-7589
<http://www.dnamed.co.kr>
E-mail: dnamed@hanmail.net

제조원

